
TIME WARP



Quick Start Guide

US and foreign Patent(s) Pending

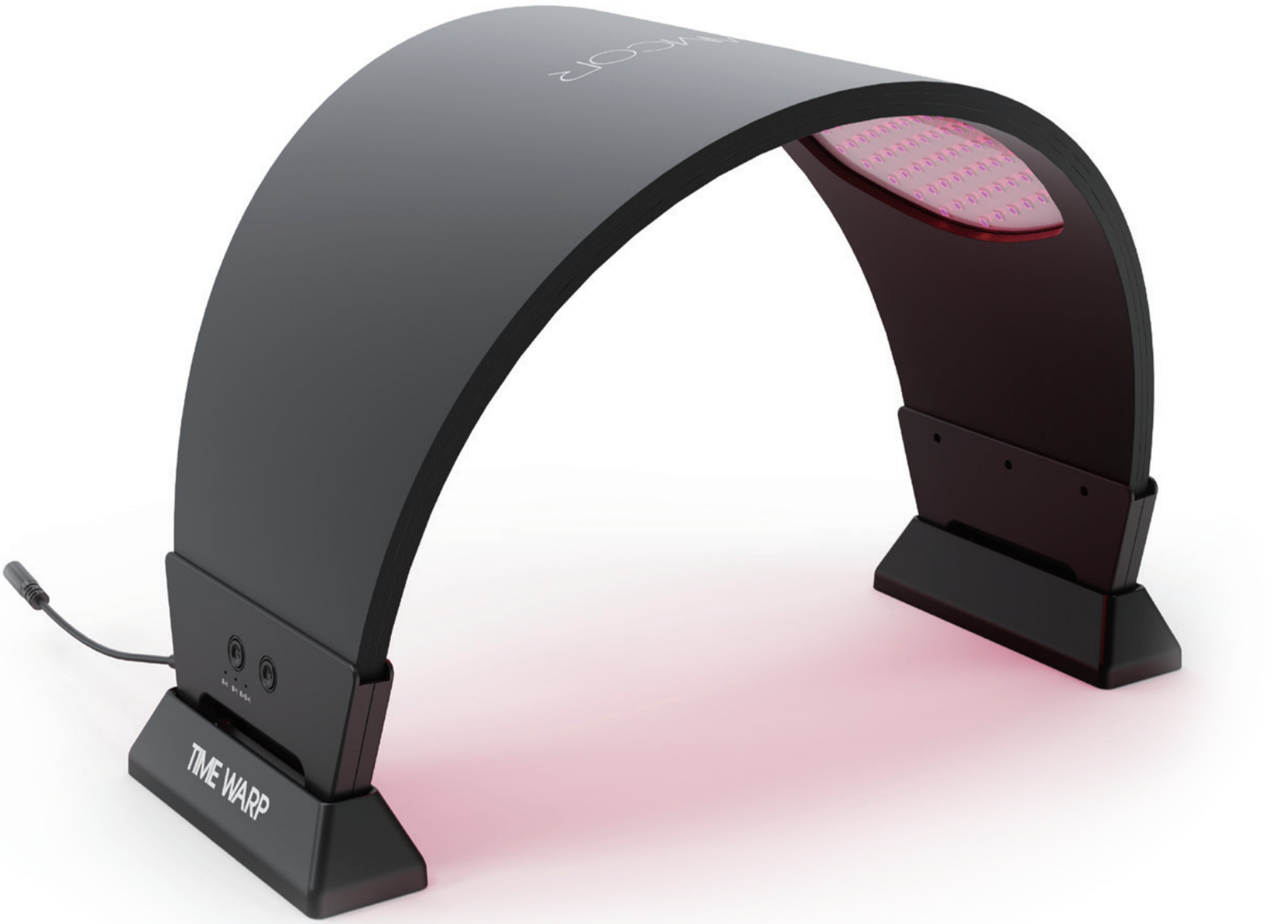

Thank you, and congratulations!

YOU JUST MADE THE BEST INVESTMENT INTO YOUR CAREER

Your TIME WARP is the perfect investment for those who seek to provide top-tier skin care treatments. This advanced therapy light, featuring a Glamcor LED dome panel, is a versatile and powerful addition to any professional skincare practice. Take a moment to review the simple instructions provided here, and you'll be ready to elevate your work with precision and style.

THANK YOU

We cannot thank you enough for choosing GLAMCOR to be part of your profession. We work really hard to constantly improve our products, and you can always expect something new from us. It is our commitment to you, and we are honored to be part of your professional career.

STAY INTOUCH

Share your success with us! You can reach us here:
Email: realperson@glamcor.com
Website: www.glamcor.com
Instagram: @glamcorofficial
Facebook: www.facebook.com/GLAMCOR

NEED SERVICE OR SUPPORT? GREAT!

We have legendary service supported by a team that operates almost around the clock. No matter where you are, send us an email sowe can log your support issue in our system. Include your phone number so we can text or call you if you would like!
Email: realperson@glamcor.com

01 What's inside?

Lightweight US adapter with 2 meter cable





Black goggles



Instructions for use

Stands for each side of mask



Stands for each side of mask

TIME WARP LED therapy mask

Product Description


Introducing the TIME WARP LED Panel Therapy Light, a professional-grade skincare tool engineered to meet the demanding needs of beauty and wellness experts. This advanced therapy light, featuring a Glamcor LED dome panel, is a versatile and powerful addition toany professional skincare practice.

02

Customizable Work Modes for Tailored Treatments


Red + Infrared Mod:

Ideal for anti-aging and deep tissue rejuvenation therapies.

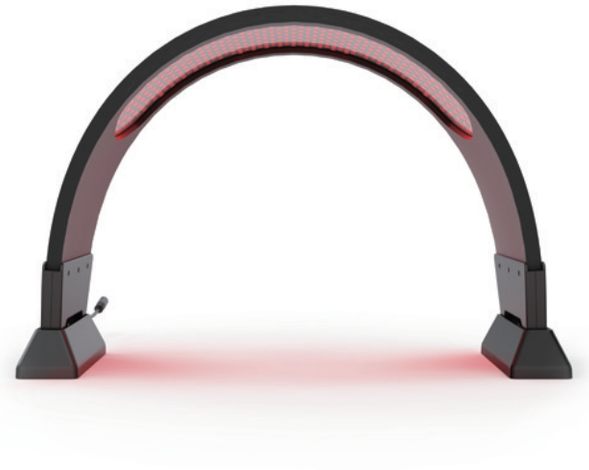

Blue + Infrared Mode:


Optimal for comprehensive acne treatment and skin calming.

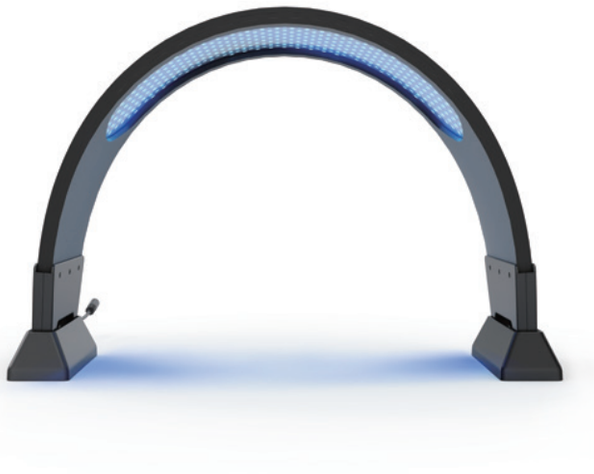

Tri-Light Therapy Mode:

Combines Red, Blue, and Infrared lights, offering a full-spectrum solution for various skin concerns



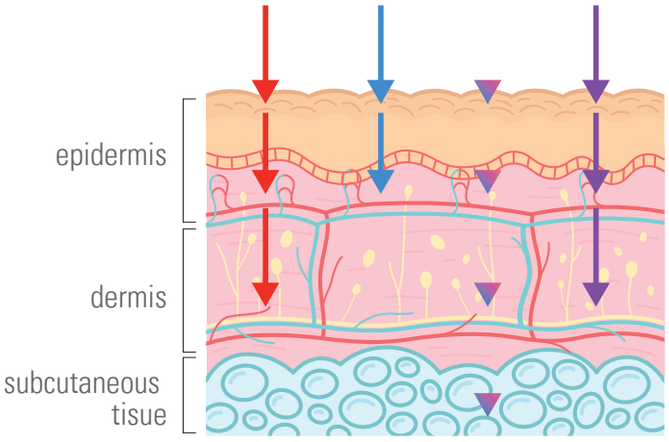

03 Time Wrap Specifications

Produce name	Time Warp
Dimensions	785mm*283mm*19mm
Weight	
LED Bulbs	430 pcs five-cores LEDs
LED Wavelengths	465nm,630nm,830nm
Output Power Density	Red 630nm: 15mw/cm2;
	Blue 465nm: 15mw/cm2;
	infrared 830nm: 15mw/cm2.
Power Supply	DC port
Input	110-260V
Output	12V=4A
Operating Temperature	-10°Cto 35C/14Fto 95F
Storing/Transport Temperature	20°C to60C/68Fto 140F
Operating Humidity	30% to 80% without condensation
Storage/Transport Humidity	30% to 90% without condensation

04 Description for Buttons



❶ Power ON/OFF

❷ Three Modes:

1.R+I: Red + Infrared Mode

2.B+I: Blue + Infrared Mode

3.R+B+I: Tri-Light Therapy Mode

05 How to use?



❶ Prepare the stand on each side



❷ Insert the panel into the stand



❸ Plug in the power



❹ Select your treatment of either R+I,B+I or R+B+I light therapy



❺ Your 10-minute treatment will begin automatically



❻ Turn off the panel by pressing the power button

Recommendation of Use

Please read the complete instruction manual before using the Glamcor LED dome panel. It will give you a better understanding of how the product works, and you'll be able to enjoy the best treatment results.

Warning:

1. This device is not for use over the eyes and is not for use on eyelids.Do not look into the light-emitting diodes when using this device.
2. Do not attempt to operate this equipment if there is evidence showing that this product is damaged, or you suspect that this product is damaged.Ask for technical support if you find any damage to the product.
3. Only authorized personnel may service this equipment, contact GLAMCOR for servicing or repairing.
4. Do not put body weight on the device. The device cannot withstand a human weight or other heavy external pressure. Do not place the device under clothing, blankets, sheets or towels.
5. Please connect the DC-port to the adaptor before plugging the adapter into to the power supply.
6. Please do not place the Glamcor LED dome panel where it may fall or be pulled into water or other liquid. If the device falls into water or liquid, please unplug immediately.

07 Recommendation of Use

7. Do not use the Glamcor LED dome panel while you are bathing or taking a shower.
8. Keep the device away from heat and hot surfaces.

9.Do not attempt to disassemble the Glamcor LED dome panel, which may cause damage to the device and void the warranty.
10. Do not let children under 18 years of age use the Glamcor LED dome panel without close supervision.

11.If the device overheats, turn it off immediately.
12. Keep the device away from flammable liquids and aerosols.
13. Use of accessories other than those specified or provided by the manufacturer of this equipment could result in damage or harm to the user.
14. We do not recommend use of the device if you have had a sun burn within two weeks.
15. If you have received any facial fillers or botox within the past month please discuss with your healthcare provider before using the panel.
16.If you have had skin cancer, please discuss with your healthcare provider before using the panel.

09 Intended Use and Instructions for Use

Intended Use

The Glamcor LED dome panel is intended to emit light in the red region of the spectrum, specifically indicated to treat mild to moderate wrinkles and signs of aging. It also emits light in the blue region of the spectrum, specifically indicated to reduce the bacteria that causes acne and blemishes and improves cell regeneration, skin tone and texture.

Instructions for Use

Check for these situations before use:
Please determine your skin type according to the table below before use, as it will provide more effective instructions for your choice. The device should only be used by Fitzpatrick Skin Types –IV:

The Fitzpatrick Scale



TypeI	Type II	Type III	Type IV	Type V	Type VI
Light, Pale Whlte	Light, Fair	Medium, White to Olive	Olive, Moderate Brown	Brown, Dark Brown	Black, Very Dark Brown to Black

Note: The "X" means not suitable. Please don't use the device if you have Fitzpatrick Skin in Type V or Type VI.

11 IMaintenance

Maintenance Cleaning

The panel is reusable and needs to be cleaned after each use.
Please clean the panel after each treatment, using the following method:
1. Use a clean, damp cotton cloth (use water, not other cleaning products) to wipe down the enclosure.
2. Do not oversaturate the cloth; you do not want any liquid to enter the interior of the device.
3. Note: Check the device whether some visible contaminants still exist. If there are visible contaminants on the device, repeat the above steps.
4. Use a clean cotton cloth lightly moistened with alcohol (70%-75%) to clean the enclosure. Like before, be careful to not oversaturate the cloth to ensure that no liquid enters the interior of the device
5. Visually check the cleanliness of the device and make sure there are no dust, impurities, or any other Pollutant. If necessary, repeat the above step.
6. Drying: please dry the device in the air naturally.
7. Note: Never use abrasive cleaning agents, and never immerse the device in water or other cleaning liquids to avoid damage to the device.

Warning:

a) DO NOT immerse the product in liquid.
b) DO NOT dismantle the unit by yourself.
c) DO NOT drop the device or apply external pressure to the unit since this may damage or weaken the LEDs and other components.
d) If the cord or device is damaged or does not operate normally, contact the distributor. This device may only be serviced, repaired, and disassembled through GLAMCOR directly. All units must be returned to the manufacturer for repair, including Warranty repair and Out-Of-Warranty repair, please do not disassemble or repair the returned units by yourself.

Light Sensitivity Test

Follow the instructions below to check if you are light-sensitive:

1. First, turn on the panel (see instructions).

2. Then, apply the panel to the facial area, making sure the light is shining on the skin, and wear for 10 minutes.

3. Wait for 24 hours. If you experience pain or discomfort, or redness that lasts for more than 2 hours, you might be light-sensitive. Please consult a doctor for advice.

Note:

It is recommended to have a 10-minute treatment, once daily, 3 times a week. Over-treatment is not advised. Do not use any makeup or cosmetics that contain chemical or acidic components that may affect light output and lead to a decreased or ineffective treatment.

08 Contradictions and Adverse Reactions

Contradictions

(Do not use If any of the following are applicable to you)

1. Do not use the device if you have the Fitzpatrick Skin in Type V and Type VI.
2. Do not use this device if you are photosensitive.
3. Do not use on infants or children.
4. Do not use if you have a medical history of seizures that are triggered by light.
5. Please do not use Glamcor LED dome panel on the thyroid gland in the neck area.
6. Do not use if you have severe acne.
7. This device should not be used by people who are:

-Under the age of 18

- Pregnant or breastfeeding
- Using photosensitive drugs

Adverse Reactions

If you have any adverse reaction from using this device, stop using it immediately.

Note: If you experience any discomfort, such as skin discoloration, skin damage, or any other concerns develop, stop using the device immediately and consult your physician.

Safety Instructions

Proper Use

This fan is designed for use in dry indoor rooms. It is not suitable for use in damp rooms, e.g. bathrooms. The fan is intended for air circulation. Only use this product as described in this manual.

DO NOT SUBSTITUTE A POWER SUPPLY OTHER THAN THE INCLUDED APPROVED GLAMCOR BRANDED POWER SUPPLY OR GENUINE REPLACEMENT.

Risk Of Injury For Children, Person, And Animals

- Keep children away from packaging - risk of suffocation!
- Keep unsupervised children away from the fan.
- Do not insert objects into the rotating parts of the fan.
- Do not touch the blades.
- Make sure that the fan is unplugged from the power source before performing any action that requires touching the fan blades, such as inspections or cleaning.

Location

- Place the fan at least 8” or 20 cm from flammable materials.
- Always place the fan on an even and stable surface to prevent tipping over.
- Do not subject the device to extreme temperatures, severe changes in temperatures, direct sunlight, or moisture.
- Do not allow the device to fall, and do not subject it to strong impacts.
- Keep the fan well away from heat sources such as stovetops or ovens.
- Ensure that the device is not able to fall in water and cannot get wet.

Safety Instructions (con't)

Electricity - Risk Of Electric Shock

- Keep the fan away from humidity.
- Connect the fan only to outlets corresponding with the technical data of the lamp.
- Check the fan before each use for damage.The fan is not to be used when it shows visible damage.
- Do not take the fan apart and only have it repaired at your local electrical store.
- Before cleaning or if the fan is not to be used for a long period, unplug the fan from the main power supply.
- During storms; unplug the fan from the main power supply. When removing the plug from the socket, always pull on the plug and not the cable.
- A damaged cable must only be repaired at your local electrical store.
- Do not use the fan if it is damaged or if the connecting cable shows any kind of visible damage.

Damages On Device And Materials

- Do not cover fan, for example, with drapes or anything else. This could lead to fire.

Changing Fans

- Do not attempt to change fans.

Care / Maintenance

- Keep the fan away from humidity.
- Unplug the fan from the mains power supply.
- Dust off the fan if necessary with a dry cloth.

Warranty / Guarantee

Important

Thank you for purchasing this product. It has been made to demanding high-quality standards and is guaranteed for domestic use against manufacturing faults for a period of 12 months from the date of purchase.

This guarantee does not affect your statutory rights. If your product fails due to a defect in material or workmanship during this period, please contact GLAMCOR by email (realperson@glamcor.com) or return it to the place of purchase.

Normal wear and tear is not covered under the guarantee. Any guarantee is invalid if the product has been misused or subject to neglect or an attempted repair other than by our own service center.

Due to continuous product improvement, we reserve the right to change the product specification without prior notice.

After Sales Customer Service Department:
GLAMCOR
110 Wall Street Floor 5
New York, NY 10005
EMAIL: realperson@glamcor.com

Please retain for future reference. Colors and contents may vary.

Recycling

Disposal

The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner, and make it available for the recyclable material collection service.

Dispose of the product in an environmentally friendly manner when you decide to part with it. This device should not be disposed of with household waste. Dispose of it at a recycling center for old electrical and electronic appliances.

Please recycle where facilities exist.

Check with your Local Authority or retailer for recycling advice.

FCC Notice - USA and Canada

⚠ Warning: For Indoor Use Only

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

– Reorient or relocate the receiving antenna.
– Increase the separation between the equipment and receiver.
– Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
– Consult the dealer or an experienced radio TV technician for help.

This Class B digital apparatus complies with Canadian ICES-003.
Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.

© 2024 GLAMCOR GLOBAL LLC. Reproduction in whole or in part without written permission is prohibited. All rights reserved. GLAMCOR is a trademark of GLAMCOR GLOBAL LLC.